



## Multiple Activity Cage


Cat. No. 47420

### General

The **47420 MULTIPLE ACTIVITY CAGE** package comprises:

- an **electronic unit**, Cat. 7441
- an **I.R. Beam Cage**, which consists of an Animal Cage of clear Perspex, Cat. 7433, complete with two sets of emitter/sensor arrays for horizontal and vertical activity, respectively Cat. 7435 and Cat. 7436.

This set-up can accept up to 5 additional cages, for a total of 6.

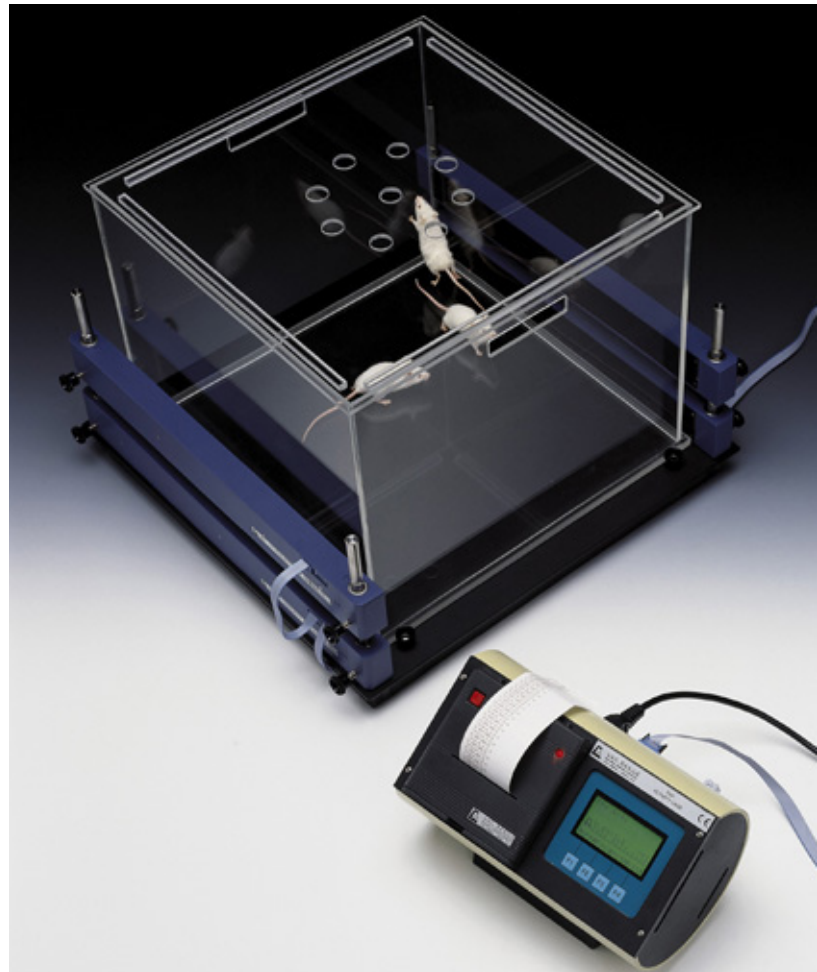

**NEW !! With dedicated software included!**

**STAND-ALONE CONTROLLER WITH EMBEDDED PRINTER FOR GLP AND DATA SAFETY**

### MAIN FEATURES

- Measures horizontal and vertical activity in rats and mice
- Psychopharmacology, in screening drugs which are potentially active on the central nervous system.
- Behavioural Sciences, in evaluating the variations of spontaneous activity after changes in environmental conditions

## Electronic Unit

The **7441** is designed to process the data originated by up to 6 **7433 cages**.

The Electronic Units incorporate a graphic display, a thermal printer and a serial port RS232 for direct connection to the PC via the software Cat. 52050 included. A serial to USB connector is also included.

The graphic display presents all available commands. The operator sets the experiment configuration via the keyboard located below the display.

The activity data are displayed at preset intervals and printed/routed to the computer according to the selected configuration. The data can be customized by adding ani-mal & experiment numbers, gender, etc.

**7441** is provided with an internal memory, capable to store the data of several experiments, to be unloaded to the PC later.

## Cage

The **7433** Cage consists of a cubicle, dimensioned 41 x 41 x 33 (h) cm, entirely made of clear Perspex, upper lid and bottom catch pan detachable for cleaning.

The cubicle rests on a sturdy base made of black Perspex, provided with four vertical notched bars of stainless steel to which the horizontal/vertical detecting systems 7435 and/or 7436 can be fastened.

The **7435** consists of two facing blocks containing an I.R. array of emitters and, respectively, sensors. The **7435** records the **horizontal activity**. A similar system, Cat. **7436**, whose height can be adjusted, assesses the **vertical activity (rearing)**.

## Data Acquisition

The electronic unit is microprocessor controlled and features direct PC output. Internally-stored data can be routed via a 9-pin D-type connector to the PC serial port (RS232).

Data output is managed by **52050-04** Data Acquisition Software Package (Windows® based), which enables the research worker to store the data into individual files, ready to be easily managed by most statistical analysis packages available on the market.

**Ask for details!**

## Ordering Information

**47420 MULTIPLE ACTIVITY CAGE**, standard package, including 7441 Electronic Unit (for up to 6 cages) and one 7433 with 7435 & 7436 emit-ter/sensor kit, cables & manual.

<b>7441</b>	Electronic Unit,
<b>7433</b>	Animal Cage complete with I.R. Beam Array
<b>7435</b>	Set of emitter/sensor arrays for horizontal activity
<b>7436</b>	Set of emitter/sensor arrays for vertical activity
<b>37400-305</b>	Package of 10 Heat Sensitive Paper Rolls
<b>7439</b>	Instruction Manual
<b>E-WP008</b>	Mains Cord
<b>52010-320</b>	USB to serial port converter
<b>52010-322</b>	Serial cable 9 to 9 pin

Set of fuses for either 230 or 115 V operation

## Physical (Dimensions & Weight)

<b>7441</b>	27x16x19 cm, Kg 2.70
<b>7433</b>	54x50x37 cm, Kg 11.80 (incl. 7435/7436)

## Bibliography

- A. Marazioti et alia: "**Somatostatin Receptors in the Ventral Pallidum/Substantia Innominata Modulate Rat Locomotor Activity**" Psychopharmacol., 181:2, 319-326, 2005
- W. Ponti et alia: "**In vivo Model for the Evaluation of Molecules Active Towards Transmissible Spongiform Encephalopathies**" Veter. Res. Communat., 28:1, 307-310, 2004
- T. Dolezal et alia: "**Guaifenesin Enhances The Analgesic Potency of Paracetamol in Mice**" Arch. Pharmacol., 366:6, 551-554, 2002
- M. Votava et alia: "**Effects of Alprazolam and Fluoxetine on Morphine Sensitization in Mice**" Physiol. Res., 51, 417-423, 2002