




Pressure Transducer Calibrator

Cat. No. 2900

General

Our **Pressure Transducer Calibrator** has been designed for easy calibration of pressure transducers.

It basically consists of a transducer holder (transducer is not provided) and a mercury manometer, provided with a sliding scale for easy zero setting, branched via a 3-way stopcock to the pressure bottle and to the transducer.

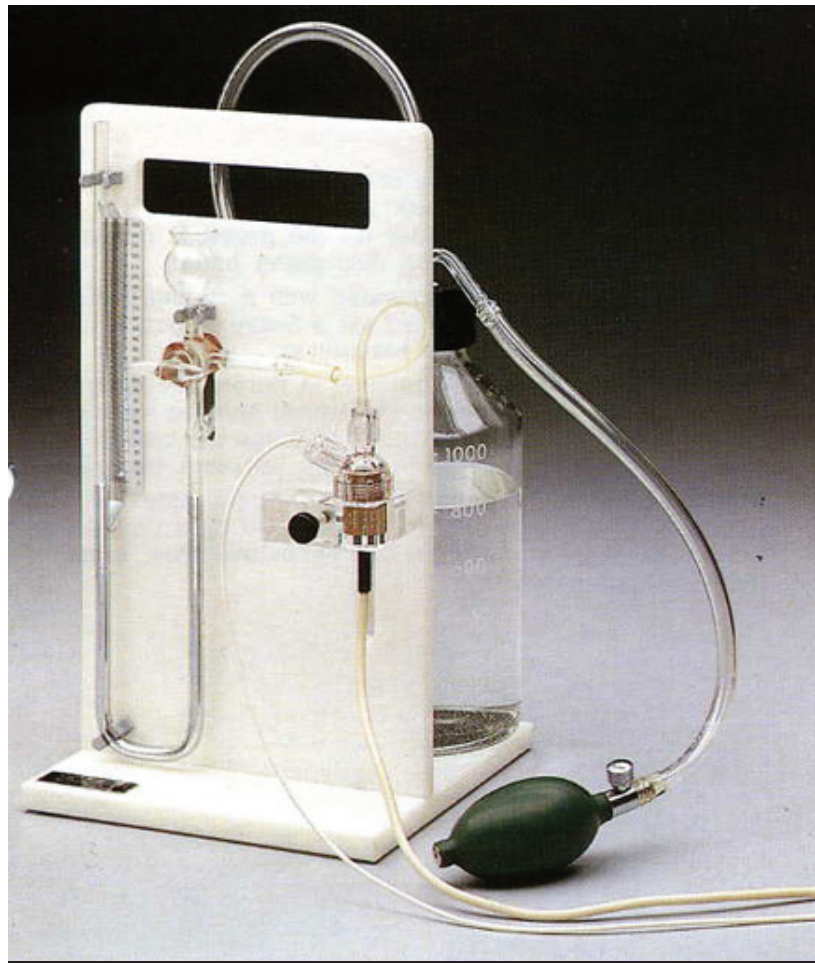

NOTE: pressure transducer is not included

Ordering Information

- 2900** Pressure Transducer Calibrator, complete
- 2901** Spare Bottle for 2900, complete with tubes
- 2902** Complete Tube for 2900
- 2903** Hand Pump

Bibliography

- J.Y. Kang et alia: "Role of Capsaicin Sensitive Nerves in Epidermal Growth Factor Effects on Gastric Mucosal Injury and Blood Flow" Gut 42: 344-350, 1998
- M. Bernardi et alia: "Spectral Analysis of Inter-cycle Heart Fluctuations in the diethyl-ether-anaesthetized or Pithed Rat Treated With Prazosin, dl-propranolol, Endothelin-1, α -r Atiopeptin and ACE-Inhibitors" J. Autonomic Pharmacol. 18 (5): 271-280, 1998